

At01 Atomic Spectroscopy & Electron Spin Resonance (Initial Report)

Matthew Evans

15th November 2018

1 Introduction

Two key atomic concepts will be investigated during this experiment: the atomic spectrum of hydrogen and the spin resonance of electrons.

Electrons occupy energy levels in atoms that orbit around the nucleus. When electrons gain energy, for example from a collision with an external electron or absorption of a quantum ‘packet’ of energy called a *photon*, they move up to the next energy level and are then said to be *excited*. After a period of time, they de-excite back to their original energy level emitting a photon equivalent to the energy difference of the two levels

$$h\nu = E_f - E_i \quad (1)$$

where $h = 6.63 \times 10^{-34}$ Js, is the Planck constant, ν is the photon frequency, E_f is the excited energy level, E_i is the original equilibrium energy level. In this experiment we will be analysing the *emitted* spectrum of hydrogen [1].

In addition, all electrons have a spin angular momentum associated with them [2]. The electron spin component is given by

$$S_z = m_s \hbar \quad (2)$$

where S_z is the electron spin angular momentum by considering the z-component, $m_s = \pm \frac{1}{2}$ refers to the spin quantum number either ‘spin up’ (positive) or ‘spin down’ (negative) respectively and $\hbar = \frac{h}{2\pi}$. This concept can be likened to an analogy of a spinning sphere [3] where the spin angular momentum is the momentum associated with the body about z-axis. The magnitude of the spin angular momentum is

$$S = \hbar \sqrt{s(s+1)} \quad (3)$$

where $s = |m_s| = \frac{1}{2}$ and is the spin quantum number. Because electrons have a charge associated with them, they will be affected as they pass through a magnetic field. This experiment will be studying this effect further and obtaining some fundamental information from this process.

2 Theory

The energy levels of a simplified diagram for the hydrogen atom is shown in Figure 1.

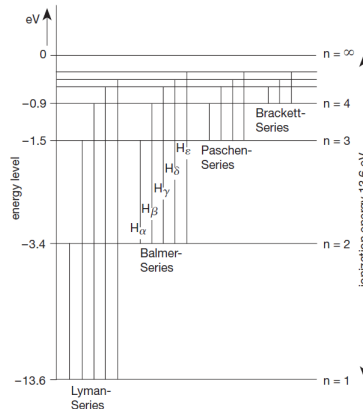


Figure 1: A simplified diagram showing the energy levels of the hydrogen atom. Where n is the principle quantum number. $n = 1$ is the ground state and $n = \infty$ corresponds to the ionisation energy required to liberate the bound electron [1].

These energy levels can be determined by using the Bohr model or by solution to the Schrödinger equation [1] and is given by

$$E_n = -\frac{Ry}{n^2} \quad (4)$$

where Ry is the Rydberg constant and has a value of 13.6 eV and n is the energy level number in the atom. This experiment will be considering transitions to the $n = 2$ energy level, i.e. the Balmer Series, for the photon energies as the first few photons in the Balmer Series lie in the visible region of light. This because the first few energy transitions for this series is in the visible region. Therefore by considering relation (1), equation (4) becomes

$$E_n = -Ry \left(\frac{1}{n^2} - \frac{1}{2^2} \right) = Ry \left(\frac{1}{2^2} - \frac{1}{n^2} \right) \quad (5)$$

For a more detailed derivation of the energies levels for the hydrogen atom using the Bohr model please refer to [4].

Electrons are negatively charged particles, therefore they have an associated magnetic moment

$$\boldsymbol{\mu} = -g \frac{e}{2m} \mathbf{S} \quad (6)$$

μ is the magnetic moment, g is the spectroscopic splitting factor, e and m are the charge and mass of the electron respectively [2]. As mentioned in Section 1 the spinning motion of the electron is effected as it passes through a magnetic field \mathbf{B} because the magnetic dipoles of the electron and field interact with each other. This means that the electron will precess about the magnetic field, \mathbf{B} , in one of two orientations as shown in Figure 2 [2].

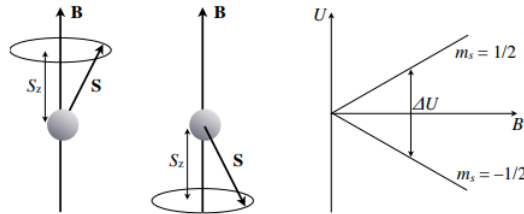


Figure 2: Left: orientation of electron ‘spin up’ and ‘spin down’ are given in a magnetic field \mathbf{B} . Right: dependence of the two energy states on \mathbf{B} is shown [2].

The magnetic energy is given by [2]

$$U = -\boldsymbol{\mu} \cdot \mathbf{B} = -\mu_z B \quad (7)$$

Where μ_z is the z-component, $\mu_z = -g \frac{e}{2m} S_z$, of the magnetic moment. By substituting (2) into (7), we now have

$$U = g \frac{e}{2m} (m_s \hbar) B = g \mu_B B \quad (8)$$

where $\mu_B = \frac{e \hbar}{2m}$ is the Bohr magneton. Therefore, from Figure 2 and (8) we can see that the difference in energy be determined by [2]

$$\Delta U = g \mu_B B \left(\frac{1}{2} - \left(-\frac{1}{2} \right) \right) = g \mu_B B. \quad (9)$$

For further information regarding this derivation arriving it equation (9), please refer to [3]. By considering the electron energy changes involved with the atomic spectroscopy mentioned earlier, we can therefore deduce that the electron can either go up or down an energy state by photon absorption or emission receptively. Hence the energy of the photon required for such a transition is given by

$$h\nu = g \mu_B B \quad (10)$$

where ν is again the frequency of the absorbed or emitted photon.

3 Method

3.1 Atomic Spectroscopy

A low pressure gas in a tube will be exposed to an electric field and this will ionise the gas releasing electrons [1]. These free electrons will collide with other electrons bound by atoms, as mentioned in Section 1. This will excite the bound electron to an excited state and will then de-excite via the emission of a photon. These emission lines will be observed through a diffraction grating as shown in Figure 3.

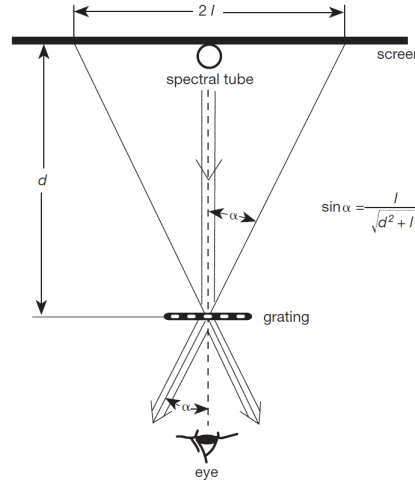


Figure 3: A plan view of the experimental set-up [1].

The first step will be to use a mercury lamp in place of the spectral tube in Figure 3 to determine the grating spacing constant d_g . The emission lines viewed through the diffraction grating correspond to the bright fringe orders, m . Hence, by using the known wavelengths of mercury given in [1] and the diffraction grating equation

$$m\lambda = d_g \sin \alpha \quad (11)$$

we can determine d_g the grating constant. Now, the mercury lamp can be removed and the spectral tube put in place as described in Figure 3. By measuring l and d the angle α can be determined along with the associated uncertainties. The emission lines viewed through the diffraction grating correspond to the bright fringe orders, m . Hence, by using the diffracting grating equation (11) we can find the wavelengths, λ of the sample for these different orders.

Once the wavelengths are determined, the energies can be obtained from the relation $E = \frac{hc}{\lambda}$. Then a graph can be plot to using (5) to determine a value for the Rydberg constant Ry .

3.2 Electron Spin Resonance

The sample for testing will be an organic compound called DPPH (1, 1-diphenyl-2-picrylhydrazil) molecule [2]. The experimental set-up is shown in Figure 4.

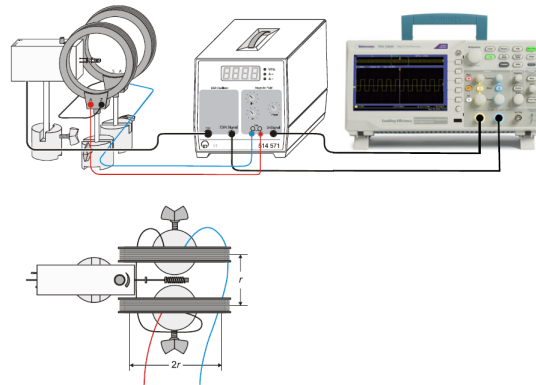


Figure 4: The experimental set-up for investigating the electron spin resonance of DPPH [2].

A series of preliminary experiments will be conducted before the main experiment commences. First, the magnetic field produced by the Helmholtz pair of coils will be measured over a range of DC currents using a Gauss meter. A relation describing the magnetic

field for this set-up in Figure 4 is given by

$$B = \mu_0 \left(\frac{4}{5} \right)^{3/2} \frac{nI}{r} \quad (12)$$

where μ_0 is the permeability of free space, n is the number of turns per coil ($n = 320$ in this case), I is the current through each coil and r is the coil radius (6.75cm) [2]. Once a range of magnetic field measurements have been obtained in conjunction with the corresponding currents, these currents can be used in (12) to obtain the magnetic field. By comparing the measured value of the magnetic field to that obtained in (12) we can justify if (12) is reliable.

The final preliminary check is to see that removing energy from a circuit will put load on the circuit and reduce the output voltage [2]. The current output from the ESR unit, see Figure 4 will be connected to a multimeter to measure the DC current. By using a passive circuit and a medium coil, the resonant frequency of the passive circuit along with that of the oscillator will be observed. When these two resonant frequencies coincide, the voltage should reach a maximum [2].

After the preliminary checks have been made, the sample can be placed in one of the three coils provided and the experiment is set-up as described in Figure 4. Once all the relevant connections have been made to the oscilloscope and the ESR unit, the frequency and current can be measured at different resonances. Once several resonance measurements for the current have been obtained, (12) can be used with the measured currents to obtain the corresponding magnetic fields. Then, these magnetic field measurements along with the corresponding frequency measurements can be used in (10) to obtain a value for g by plotting an appropriate straight-line graph or by using the individual results to obtain several values of g and then averaging to obtain a single final value.

3.3 Safety

During both experiments there are many sensitive electrical components used; such as a radio frequency oscillator, an oscilloscope and power supplies. Therefore, no food or drink should be consumed whilst working in the laboratory and any exposed metal should not be touched when any power supply is switched on [1]. In addition to this, the glass discharge tubes for the Atomic Spectroscopy experiment are fragile. So, when they require moving or changing the demonstrator or laboratory technician should be present to handle these pieces of equipment.

References

- [1] Atomic Spectroscopy worksheet.
- [2] Electron Spin Resonance worksheet.
- [3] Young, Hugh D and Freedman, Roger A, *University Physics*, 13th Edition, Chapter 41, pages 1519 - 1527.
- [4] Young and Freedman, *University Physics*, 13th Edition, Chapter 39, pages 1427 - 1437.